


ONE-CARE®

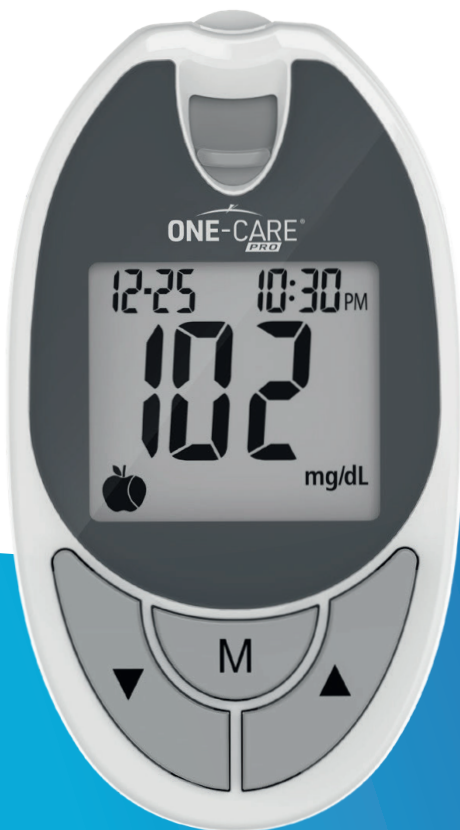

- **Blood Sampling**
- **Glucose Monitoring**
- **Hypodermic Injection Devices**
- **Personal Protective Equipment**

2022
PRODUCT
CATALOG

ONE-CARE[®] PRO

GLUCOSE MONITOR & STRIPS
FOR PROFESSIONAL USE



Auto Calibration



0.5 µL Blood Sample



Alternative Site Testing



Instant 14 Days Average




5 Second Test



500 Memory Sets



PC Download

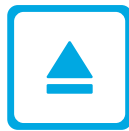

Strip Ejector

Bigger Buttons in a compact device

Easy to use




Auto Coding

Simplified way to test, easy to use, no more misreading due to incompatible code card

GOD Enzyme

0.5 µl blood sample
Highly Accurate for Capillary Whole Blood samples

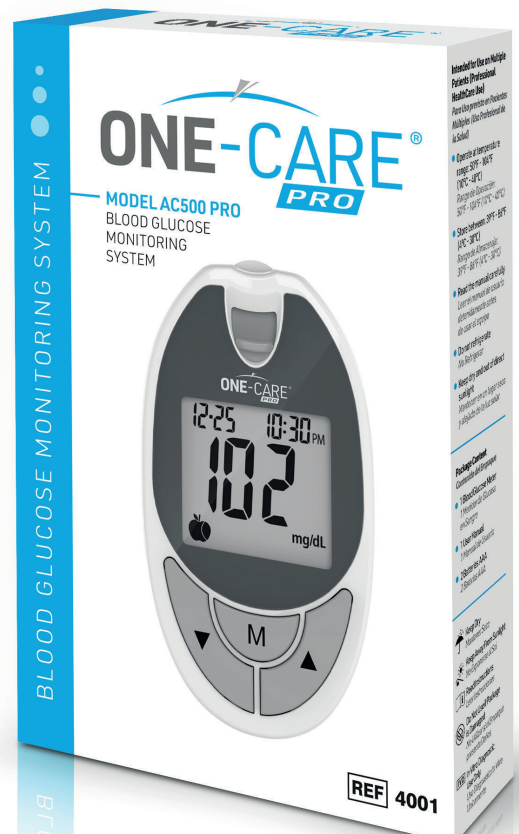

For more information direct to your smartphone

SCAN THE CODE
> GLUCOSE MONITOR & STRIPS FOR MULTI-PATIENT



Hematocrit	30-55%
Test Sample	Capillary Whole Blood
Measuring Unit	mg/dL or mmol/L
Measuring Range	20 - 600 mg/dL (1,1 - 33,3 mmol/L)
Average	14/07/30 days
Reminder Alarm	4 sets
External Output	UART with mini USB socket (Optional) Bluetooth (Optional)
Power Supply	2 AAA batteries
Battery Life	Approximately 1,000 tests
Dimensions	100 (L) x 57 (L) x 24 (A) mm
Weight	58 g without batteries

Test Strip Ejector

Reduce the possibility of cross-contamination due to bloodborne pathogens



Test Strip

Longer shelf life 24 months



PRODUCT ORDERING INFORMATION

REF.	PRODUCT	PCS per BOX	BOXES per CASE
4001	ONE-CARE® PRO GLUCOSE METER FOR PROFESSIONAL USE	1	80
4002	ONE-CARE® PRO GLUCOSE STRIPS FOR PROFESSIONAL USE	50	560
4003	ONE-CARE® PRO CONTROL SOLUTION LEVEL 1	1	900
4004	ONE-CARE® PRO CONTROL SOLUTION LEVEL 2	1	900
4005	ONE-CARE® PRO CONTROL SOLUTION LEVEL 3	1	900
4010	ONE-CARE® CONTROL SOLUTION LEVEL 1 & 3 (LOW & HIGH)	2	900

ONE-CARE

SAFETY LANCETS

OPTI



Ergonomic
Ergonomic body for comfortable grip

Easy
Wide cap for easy twist off
















Comfortable
Dual spring safety activation mechanism for enhanced patient's comfort



Every product detail designed with ease of use and comfort in mind

PRODUCT ORDERING INFORMATION

Available in seven color-coded versions to meet different test blood flow requirements and address a variety of patient skin types.

REF.	VERSION	DIAMETER (G) / DEPTH (mm/in)	COLOR	PACKAGING
1101	Micro Flow	30G / 1.5 mm 		100 pcs / box 2,000 pcs / case
1102	Mini Flow	28G / 1.8 mm 		100 pcs / box 2,000 pcs / case
1103	Low Flow	26G / 1.8 mm 		100 pcs / box 2,000 pcs / case
1104	Normal Flow	23G / 1.8 mm 		100 pcs / box 2,000 pcs / case
1105	Universal Flow	21G / 1.8 mm		100 pcs / box 2,000 pcs / case
1106	High Flow	21G / 2.4 mm 		100 pcs / box 2,000 pcs / case
1107	Max Flow <small>Blade</small>	18G / 1.8 mm 		100 pcs / box 2,000 pcs / case



For more information direct to your smartphone

SCAN THE CODE
> SAFETY LANCETS

ONE-CARE[®] PRO

PUSH BUTTON SAFETY LANCETS

Ergonomic


Ergonomic T-Shape body for convenient grip

Comfortable

Dual spring safety activation mechanism for enhanced patient's comfort





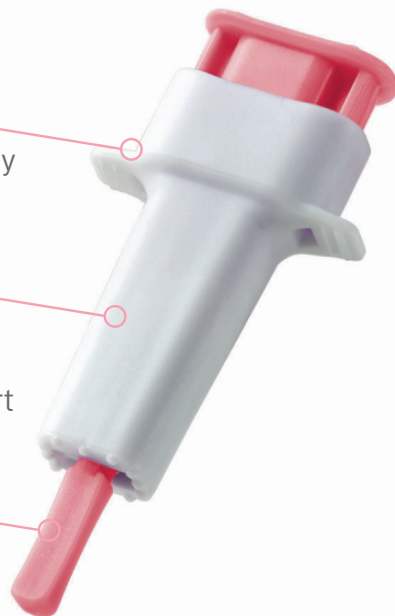







Easy

Long cap for easy twist off and pull out



PRODUCT ORDERING INFORMATION

Available in six color-coded versions to meet different test blood flow requirements and address a variety of patient skin types.

REF.	VERSION	DIAMETER (G) / DEPTH (mm/in)	COLOR	PACKAGING
1011	Micro Flow	30G / 1.2 mm 		100 pcs / box 2,000 pcs / case
1012	Mini Flow	28G / 1.5 mm 		100 pcs / box 2,000 pcs / case
1013	Low Flow	25G / 1.6 mm 		100 pcs / box 2,000 pcs / case
1014	Normal Flow	21G / 1.8 mm 		100 pcs / box 2,000 pcs / case
1015	High Flow	21G / 2.4 mm 		100 pcs / box 2,000 pcs / case
1016	Max Flow	21G / 2.8 mm 		100 pcs / box 2,000 pcs / case



For more information direct to your smartphone

SCAN THE CODE
> PRO SAFETY LANCETS

ONE-CARE[®]

PLUS

CONTACT-ACTIVATED SAFETY LANCETS

Ergonomic
Ergonomic body for comfortable grip

Reliable
Reliable safety activation mechanism

Comfortable
Wide cap for easy twist off



PRODUCT ORDERING INFORMATION

Available in six color-coded versions to meet different test blood flow requirements and address a variety of patient skin types.

REF.	VERSION	DIAMETER (G) / DEPTH (mm/in)	COLOR	PACKAGING
1001	Micro Flow	30G / 1.5 mm		100 pcs / box 2,000 pcs / case
1002	Mini Flow	28G / 1.8 mm		100 pcs / box 2,000 pcs / case
1003	Low Flow	25G / 1.8 mm		100 pcs / box 2,000 pcs / case
1004	Normal Flow	21G / 1.8 mm		100 pcs / box 2,000 pcs / case
1005	High Flow	21G / 2.4 mm		100 pcs / box 2,000 pcs / case
1006	Max Flow	21G / 2.8 mm		100 pcs / box 2,000 pcs / case



For more information direct to your smartphone

SCAN THE CODE
➤ PLUS SAFETY LANCETS

ONE-CARE[®]

SAFETY PEN NEEDLES



Safety

Clear safety shield automatically locks after use to protect the needle

Simple


Color indicator becomes invisible after engagement of safety mechanism

Ergonomic Design

Ergonomic main housing with ribs for comfortable grip



Universal

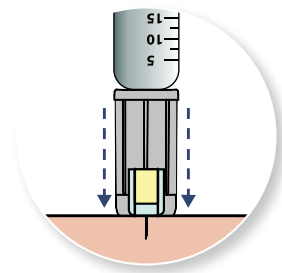Compatible fit with most pen injectors



PRODUCT ORDERING INFORMATION

Available in two color-coded product versions for the right fit for every patient and skin type.

REF.	VERSION	DIAMETER (G) / DEPTH (mm/in)	COLOR	PACKAGING
3101	Micro	30G x 5mm (3/16 in)		100 pcs / box 2,400 pcs / case
3102	Normal	30G x 8mm (5/16 in)		100 pcs / box 2,400 pcs / case




For more information direct to your smartphone

SCAN THE CODE
> SAFETY PEN NEEDLES

ONE-CARE®

INSULIN SAFETY SYRINGES

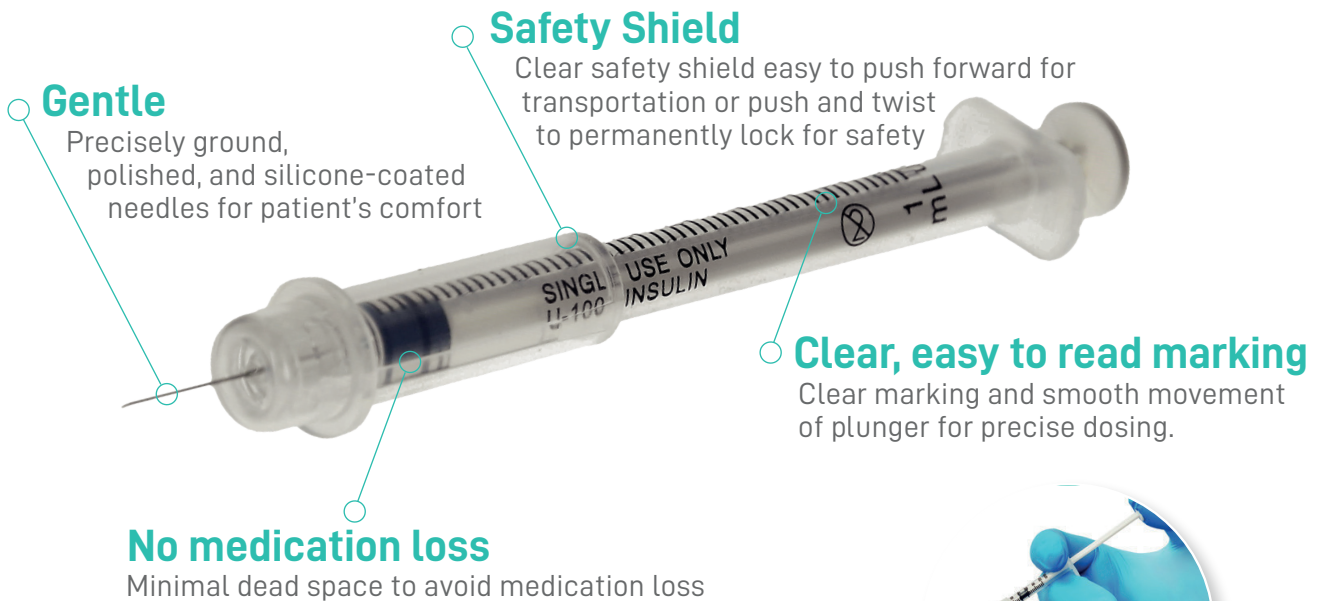

PRODUCT ORDERING INFORMATION

Available in two versions for the right fit for every patient and dose requirement.

REF.	SIZE	NEEDLE (DIAMETER x LENGTH)	COLOR	PACKAGING
2101	1 mL	29G x 1/2" 0.3mm x 12.7mm	■	100 pcs / box 500 pcs / case
2102	0.5 mL	29G x 1/2" 0.3mm x 12.7mm	■	100 pcs / box 500 pcs / case



For more information direct to your smartphone

SCAN THE CODE
> INSULIN SAFETY SYRINGES

Safety Syringes and Safety Needles

Luer-Lock Syringes 1ml with Safety Needles 23G x 1" (25mm)



REF 2110

1 ML 23G x 1" (25mm)

Easy, one-hand activation of safety shield

Audible click indicates safety mechanism activation

Locking mechanism ensures needle is secure and completely locked in the safety shield

PRODUCT/APPLICATION	PRODUCT DESCRIPTION	REF	BOX / CASE
Syringes with Safety Needles for Injections	ONE-CARE® Luer-Lock Syringes 1ml with Safety Needles 23G x 1" (25mm)	2110	100 pcs / bx 12 bx / cs

Luer-Lock Syringes 1ml with Safety Needles 25G x 1" (25mm)



REF 2112

1 ML 25G x 1" (25mm)

Easy, one-hand activation of safety shield

Audible click indicates safety mechanism activation

Locking mechanism ensures needle is secure and completely locked in the safety shield

PRODUCT/APPLICATION	PRODUCT DESCRIPTION	REF	BOX / CASE
Syringes with Safety Needles for Injections	ONE-CARE® Luer-Lock Syringes 1ml with Safety Needles 25G x 1" (25mm)	2112	100 pcs / bx 12 bx / cs

ONE-CARE®

Safety Syringes and Safety Needles

SAFETY NEEDLES

Safety Needles 23Gx1" (25mm)

REF 2111



Reliable safety mechanism integrated with the needle

Convenient, one-hand use

Compatible with all Luer-Lock and Luer-Slip Syringes as per ISO 80369

Bevel oriented needle for easy injection

Compliant with performance requirements of ISO 7864

Color coding according to ISO 6009

PRODUCT/APPLICATION	PRODUCT DESCRIPTION	REF	BOX / CASE
Safety Needles for Injections	ONE-CARE® Safety Needles 23G x 1" (25mm)	2111	100 pcs / bx 20 bx / cs

Safety Needles 25Gx1" (25mm)

REF 2113



Reliable safety mechanism integrated with the needle

Convenient, one-hand use

Compatible with all Luer-Lock and Luer-Slip Syringes as per ISO 80369

Bevel oriented needle for easy injection

Compliant with performance requirements of ISO 7864

Color coding according to ISO 6009

PRODUCT/APPLICATION	PRODUCT DESCRIPTION	REF	BOX / CASE
Safety Needles for Injections	ONE-CARE® Safety Needles 25G x 1" (25mm)	2113	100 pcs / bx 20 bx / cs

ONE-CARE®

Syringes and Needles

LUER-SLIP DISPOSABLE SYRINGES, 3-PART

REF 2001

1 ml




Specially designed 1 ml syringes for accurate dose administration

Silicone lubrication

Ensures smooth plunger movement

Safe plunger backstop

Prevents accidental withdrawal of plunger

Large and stable plunger grip

Ensures stability and comfort during injection

High-transparency barrel

Makes it easy to see syringe content

PRODUCT/APPLICATION	PRODUCT DESCRIPTION	REF	BOX / CASE
Disposable Syringes for Vaccine Injections	ONE-CARE® Luer-Slip Syringes 1ml, 3-part	2001	100 pcs / bx 24 bx / cs

LUER-LOCK DISPOSABLE SYRINGES, 3-PART

REF 2002

1 ml




Luer-Lock 1 ml syringes for secure connection with hypodermic needles

Silicone lubrication

Ensures smooth plunger movement

High-transparency barrel

Makes it easy to see syringe content

Safe plunger backstop

Prevents accidental withdrawal of plunger

Large and stable plunger grip

Ensures stability and comfort during injection

PRODUCT/APPLICATION	PRODUCT DESCRIPTION	REF	BOX / CASE
Disposable Syringes for Vaccine Injections	ONE-CARE® Luer-Lock Syringes 1ml, 3-part	2002	100 pcs / bx 24 bx / cs

ONE-CARE®

Syringes and Needles

DISPOSABLE NEEDLES

ONE-CARE® DISPOSABLE NEEDLES 23G X 1" (25mm) **REF 23100-CM**

ONE-CARE® DISPOSABLE NEEDLES 23G X 1.5" (38mm) **REF 23150-CM**



ONE-CARE® hypodermic needles for vaccine administration have precise tri-bevel design, are electro-polished and lubricated for gentle and virtually pain-free injections.

ONE-CARE® DISPOSABLE NEEDLES 21G X 1.5" (38mm) **REF 21150-CM**











ONE-CARE® Canula Needles 21G x 1.5" (38mm) are recommended for vaccine preparation as may be required by vaccine manufacturer's instructions for use.

PRODUCT/APPLICATION	PRODUCT DESCRIPTION	REF	BOX / CASE
Disposable needles for vaccine dilution	ONE-CARE® Hypodermic Needles 21G x 1.5" (38mm)	21150-CM	100 pcs / bx 100 bx / cs
	ONE-CARE® Hypodermic Needles 23G x 1.5" (38mm)	23150-CM	
	ONE-CARE® Hypodermic Needles 23G x 1" (25mm)	23100-CM	
Hypodermic needles for vaccine injections	ONE-CARE® Hypodermic Needles 23G x 5/8" (16mm)	23058-CM	100 pcs / bx 100 bx / cs
	ONE-CARE® Hypodermic Needles 25G x 1.5" (38mm)	25150-CM	
	ONE-CARE® Hypodermic Needles 25G x 1" (25mm)	25100-CM	
	ONE-CARE® Hypodermic Needles 25G x 5/8" (16mm)	25058-CM	

ONE-CARE®

GLUCOSE MONITORING SYSTEM
FOR SELF-TESTING



- 
FREE
 Auto Calibration
- 
0.5µl
 0.5 µL Blood Sample
- 
AST
 Alternative Site Testing
- 
5 sec
 5 Second Test
- 
14
 Instant 14 Days Average
- 
500
 500 Memory Sets
- 
PC
 PC Download
- 
Strip
 Strip Ejector

TEST KIT

- > 1 Blood Glucose Meter with batteries
- > 10 Blood Glucose Test Strips
- > 1 Adjustable Depth Lancing Device
- > 10 Sterilized Lancets
- > 1 Protective Carrying Case
- > User Manual, Package Inserts

Note: The contents of a Complete Kit may vary and are subject to change according to the actual purchase.



PRODUCT ORDERING INFORMATION

REF.	VERSION	PCS per BOX	BOXES per CASE
4006	ONE-CARE® GLUCOSE MONITORING KIT FOR SELF-TESTING	1	80
4007	ONE-CARE® GLUCOSE STRIPS FOR SELF-TESTING	50	560



For more information direct to your smartphone

SCAN THE CODE
> GLUCOSE METER

ONE-CARE®

Lancing Device & Lancets



Easy
Easy Depth Setting
with clear marking


"Ready to Activate"
Visual Indicator

Compact Size
Compact Size for discrete
blood sampling

Convenient
Convenient Loading System

Universal
Universal Design for
fit with most common lancets



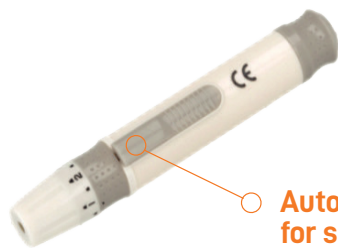





**Ultra-thin, electro-polished
lancets for comfortable sampling**



**Automated Ejector
for safe lancet removal**



PRODUCT ORDERING INFORMATION

REF.	DIAMETER (G) / DEPTH (mm/in)	DESCRIPTION	COLOR	PACKAGING
1200	10 - Depth Setting	Universal Lancing Device		1 pc/ box 100 pcs / case
1201	23G	23G Universal Blood Lancets		100 pcs / box 10.000 pcs / case
1202	26G	26G Universal Blood Lancets		100 pcs / box 10.000 pcs / case
1203	28G	28G Universal Blood Lancets		100 pcs / box 10.000 pcs / case
1204	30G	30G Universal Blood Lancets		100 pcs / box 10.000 pcs / case
1205	31G	23G Universal Blood Lancets		100 pcs / box 10.000 pcs / case
1206	32G	32G Universal Blood Lancets		100 pcs / box 10.000 pcs / case
1207	33G	33G Universal Blood Lancets		100 pcs / box 10.000 pcs / case



For more information direct
to your smartphone

SCAN THE CODE
> LANCING DEVICE & LANCETS

ONE-CARE[®]

PEN NEEDLES



Gentle

Precisely grinded three bevel tip, electro-polished and lubricated canulas for enhanced comfort

Optimal Insulin Flow

Thin- and Extra-Thin-Wall Technology

Universal Fit

ISO 11608-2 compliant and universal thread design for compatibility with all type A insulin pen injectors

PRODUCT ORDERING INFORMATION

Available in four color-coded versions to meet different test blood flow requirements and address a variety of patient skin types.



REF	DIAMETER (G) x DEPTH (mm/in)	COLOR	PACKAGING
3011	33G x 4mm	Pink	100 pcs / box 5,000 pcs / case
3012	32G x 4mm	Light Green	100 pcs / box 5,000 pcs / case
3013	32G x 5mm	Brown	100 pcs / box 5,000 pcs / case
3014	32G x 6mm	Yellow	100 pcs / box 5,000 pcs / case
3015	32G x 8mm	Grey	100 pcs / box 5,000 pcs / case
3016	31G x 4mm	Light Blue	100 pcs / box 5,000 pcs / case
3017	31G x 5mm	Purple	100 pcs / box 5,000 pcs / case
3018	31G x 6mm	Dark Blue	100 pcs / box 5,000 pcs / case
3019	31G x 8mm	Orange	100 pcs / box 5,000 pcs / case



For more information direct to your smartphone

SCAN THE CODE
> PEN NEEDLES

ONE-CARE[®]

ALCOHOL PREP PADS



ONE-CARE™ Alcohol Prep Pads are saturated with 70% Isopropyl Alcohol for optimal antiseptic action (anti-bacteria, anti-virus). They are intended for topical use in preparation of skin prior to an injection or blood collection, cleaning of minor cuts and scrapes, or cleaning of medical equipment and supplies. Reliable 4-layer packaging provides air-tight seal and prevents leakage or dry-out.

FEATURES

- Soft, absorbent, nonwoven pads saturated with 70% v/v Isopropyl Alcohol
- For Individual, Professional and Hospital Use
- Reliable 4-layer packaging
- Compact Box of 100 pads
- Attractive Price

BENEFITS

- Exceeds CDC's recommended 60% alcohol content. Optimal antiseptic action helps to reduce bacteria.
- Intended for use in preparation of skin prior to an injection or blood collection, cleaning of minor cuts and scrapes, or cleaning of medical equipment and supplies
- Air-tight seal prevents leakage or dry-out
- Safe to use in medical setting
- Lower cost of warehousing and shipping
Financial value / cost savings for quality and reliable product as compared to major national brands

PRODUCT ORDERING INFORMATION

CODE	DESCRIPTION 1	DESCRIPTION 2	SIZE	PACKAGING
1020	ONE-CARE™ Alcohol Prep Pads	Saturated with 70% Isopropyl Alcohol	M	Box 100 pcs Case 50 Boxes (5,000 pcs) Pallet 6 cases (6,300 boxes, 630,000 pcs)



For more information direct to your smartphone

SCAN THE CODE
> ALCOHOL PREP PADS

Nitrile Examination Gloves



MEDICAL GRADE, POWDER-FREE, TEXTURED EXAMINATION GLOVES

Product Description

Type of the glove

Non-sterile, powder free, examination and protective glove for single use

Material

Nitrile

Colour

Blue

Shape

Ambidextrous

Cuff

Beaded

External surface

Finger textured 35 mm. From Tip

Internal surface

Polymerized + chlorinated

Packaging

10 x 100 pcs



Physical Properties

DIMENSIONS

WIDTH (in)

80±10

(6-7)

LENGTH (in)

≥9.44

≥240

minimum

THICKNESS (mm)/(mil)

(Single wall)

minimum

ELONGATION AT BREAK (%)

minimum

FORCE AT BREAK (N)

minimum

	S	M	L	XL
WIDTH (in)	3.14 ± 0.39	3.74 ± 0.39	4.33 ± 0.39	>4.33 ± 0.39
WIDTH (mm)	80 ± 10	95 ± 10	110 ± 10	>110 ± 10
Tolerances	(6-7)	(7-8)	(8-9)	(9-10)
LENGTH (in)	≥ 9.44	≥ 9.44	≥ 9.44	≥ 9.44
LENGTH (mm)	≥ 240	≥ 240	≥ 240	≥ 240
THICKNESS (mm)/(mil)	M - Finger: min 3.15 mil min 0.08 mm		M - Palm: min 2.36 mil min 0.06 mm	
(Single wall)	M - Cuff: min 1.97 mil min 0.05 mm			
THICKNESS (mm)/(mil)	M - 500 (ISO 11193-1)		M - 400 (ISO 11193-1)	
ELONGATION AT BREAK (%)	M - 500 (ISO 11193-1)		M - 400 (ISO 11193-1)	
FORCE AT BREAK (N)	6.0 (EN 455-2) / 7.0 (ISO 11193-1)			

Manufacturing and Safety standards


AQL	Manufacturing final release: G-I inspection level AQL 1.5 in accordance with ISO 2859-1
US FDA Registration	Class 1, Product Code LZA, 510(k) K172930, device listing D320597
CE classification	Class I - Medical Device Regulation (EU) 2017/745 Category III - Personal Protective Equipment (Regulation) (EU)
Compliances	EN 455-1, EN 455-2, EN 455-3, EN 455-4 EN ISO 15223-1 EN 1041 ISO 11193-1 EN ISO 13485 EN ISO 374-1 (Type B), EN 374-2, EN 374-4, EN ISO 374-5 EN 16523-1 EN 420
Viral test	Test in accordance with ASTM F1671 & ISO 16604
Cytostatics permeation	Test in accordance with ASTM D6978
Chemical substances permeation	Test in accordance with EN 16523-1
Food contact	Declaration of conformity for food contact in accordance with Regulation (EC) No 1935/2004 and Overall Migration Test in accordance with European Commission Directive 93/11/EEC and Council of Europe Resolution AP (2004) 4
Shelf life	3 years
Storage Instruction	Keep out of direct sunlight. Store in a cool, dry place in temperature 5-35° C. Keep away from sources of ozone and ignition

PRODUCT/APPLICATION	PRODUCT DESCRIPTION	REF	BOX / CASE
Personal Protection	NITRYLEX Nitrile Examination Gloves, Powder-Free, Blue	S: RD30019902 M: RD30019903 L: RD30019904 XL: RD30019905	100 pcs / bx 10 bx/ case

For more details, contact our sales department


Particulate Respirators Cone-Design Type N95

N95 CONE-DESIGN DISPOSABLE RESPIRATOR, NIOSH CERTIFIED

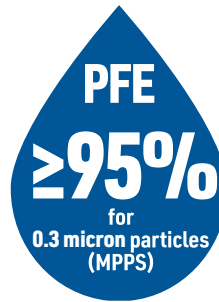

- Contour Fit Design
- Adjustable Nosepiece
- Comfortable Nose Foam
- Easy Breathing Speaking
- NIOSH N95 APPROVED

MS6115L Cone-Design Type N95

Mask Type	Standards	Filtration Efficiency*
Respirator Mask (PFE Standard)		N95
	USA: NIOSH (42 CFR 84)	PFE ≥ 95% for 0.3-micron particles

* (PFE - Particle Filtration or BFE - Bacteria Filtration Efficiency)

Respiratory Mask of Class N95 is suitable for filtering of the most penetrating particles of 0.3-micron size that can be found in the air.

0.1 micron: Particle Filtration Efficiency standard (PFE).

0.3 micron: Used to represent the most-penetrating particle size (MPPS), which is the most difficult size particle to capture.

3.0 microns: Bacteria Filtration Efficiency standard (BFE) for masks per EN 14683.


Name	MS6115L Cup-type N95
Colour	White
Item Number	6115L
Weight/unit (oz) (g)	±0.3oz 8.68g
Outer layer	PP Spunbond
Marking	As per NIOSH
Soft foam nosepiece	Black PU foam
Dimensions (in)	4.72*5.51*2.09 in (±0.20in)
Dimensions (mm)	120*140*53 mm (± 5mm)
Description	N95 Particulate Respirator
Material	3-ply, PP spunbond outer layer, PP meltbown high filtration layer, PP needle-punched nonwoven inner layer
Elastic head band	Latex-free elastic-width 7mm
Inner layer	PP needle-punched nonwoven
Filter media	PP meltbown high filtration layer
Adjustable nose clip	Aluminium
Characteristics	2 CFR standard. TC#84A-5530 95% filtration effective against particulate aerosols, free of oils, Aluminum nose-clip, Latex-Free elastic ultrasonic welded to cup

		Boxes	Cases
ORDERING INFORMATION	Materials	350 gsm white cardboard	Brown 5-ply single wall
	Dimensions (in) (mm)	5.31*4.92*7.87 (in) 135*125*200 (mm)	25.59*11.42*16.53 (in) 650*290*420 (mm)
	6115L	Gross Weight: (lb) 12.46lb/cs	(kg) 5.65kg/cs

Package 20 masks/polybag/box * 20 boxes = 400 masks / case

Particulate Respirators Foldable-Design Type N95

N95 FOLDABLE-DESIGN RESPIRATOR, NIOSH CERTIFIED

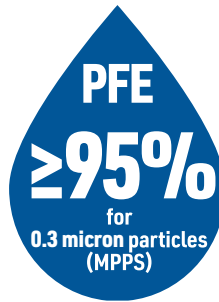

- Soft smooth lining
- Considerable fit design
- adjustable nosepiece
- Comfortable nose foam
- Easy breathing speaking

Mask Type	Standards	Filtration Efficiency *
Respirator Mask (PFE Standard)	USA: NIOSH (42 CFR 84)	N95 PFE ≥ 95% for 0.3-micron particles

* (PFE - Particle Filtration or BFE - Bacteria Filtration Efficiency)

Respiratory Mask of Class N95 is suitable for filtering of the most penetrating particles of 0.3-micron size that can be found in the air.

0.1 micron: Particle Filtration Efficiency standard (PFE).

0.3 micron: Used to represent the most-penetrating particle size (MPPS), which is the most difficult size particle to capture.

3.0 microns: Bacteria Filtration Efficiency standard (BFE) for masks per EN 14683.

MS8225 Foldable-Design N95


Name	MS8225 Folded N95
Colour	White
Item Number	8225
Weight/unit (oz) (g)	±0.34oz 9.6g
Outer layer	PP needle-punched nonwoven
Marking	As per NIOSH
Soft foam nosepiece	Black rectangular foam
Dimensions (in)	4.53*6.30 in (± 0.20in)
Dimensions (mm)	115*160 mm (± 5mm)
Description	N95 Particulate Respirator
Material	4-plys, PP spunbond outer layer, PP meltbown high filtration layer, PP needle-punched nonwoven inner layer Inner layer 19g, PP Spunbond layer
Elastic head band	Spandex headstraps
Inner layer	PP needle-punched nonwoven
Filter media	PP meltbown high filtration layer
Adjustable nose clip	Iron wire + plastic nose clip
Characteristics	Conforms to NIOSH 42 CFR standard. TC#84A-7447 95% filtration effective against particulate aerosols free of oils. Adjustable plastic nose clip, Latex-Free elastic ultrasonic welded to cup

		Boxes	Cases
ORDERING INFORMATION	Materials	350 gsm white cardboard	Brown 5-ply single wall
	Dimensions (in) (mm)	5.12*4.72*6.30 (in) 130*120*160 (mm)	24.41*11.02*13.38 (in) 620*280*340 (mm)
	8225	Gross Weight: (lb) 11.90lb/cs	(kg) 5.4kg/cs
Package 20 masks/polybag/box * 20 boxes = 400 masks / case			

ONE-CARE®

Particulate Respirators Foldable Exhalation Valve N95 Particulate Respirator

N95 FOLDABLE-DESIGN RESPIRATOR, NIOSH CERTIFIED

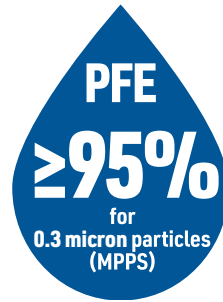

- Nonirritating and pollution free
- Easy to use & maintenance free
- Invisible / inside nose clip
- High filter material with 95% efficiency
- Adjustable headband design

MS8265 Foldable Exhalation Valve N595

Name	MS8265 Folded with valve
Colour	White
Item Number	8265
Weight/unit (oz) (g)	±0.34oz 9.6g
Outer layer	PP needle-punched nonwoven
Marking	As per NIOSH
Soft foam nosepiece	Black rectangular foam
Dimensions (in)	4.53*6.30 in (± 0.20in)
Dimensions (mm)	115*160 mm (± 5mm)
Description	N95 Particulate Respirator
Material	4-plys, PP spunbond outer layer, PP meltbown high filtration layer, PP needle-punched nonwoven inner layer
Elastic head band	Spandex headstraps
Inner layer	PP needle-punched nonwoven
Filter media	PP meltbown high filtration layer
Adjustable nose clip	Iron wire + plastic nose clip
Characteristics	Conforms to NIOSH 42 CFR standard. TC#84A-7448 95% filtration effective against particulate aerosols free of oils. Adjustable plastic nose clip, Latex-Free elastic ultrasonic welded to cup

MS8265 Folded with valve
 White
 8265
 ±0.34oz | 9.6g
 PP needle-punched nonwoven
 As per NIOSH
 Black rectangular foam
 4.53*6.30 in (± 0.20in)
 115*160 mm (± 5mm)
 N95 Particulate Respirator
 4-plys, PP spunbond outer layer, PP meltbown high filtration layer, PP needle-punched nonwoven inner layer
 Inner layer 19g, PP Spunbond layer
 Spandex headstraps
 PP needle-punched nonwoven
 PP meltbown high filtration layer
 Iron wire + plastic nose clip
 Conforms to NIOSH 42 CFR standard. TC#84A-7448 95% filtration effective against particulate aerosols free of oils. Adjustable plastic nose clip, Latex-Free elastic ultrasonic welded to cup

Mask Type	Standards	Filtration Efficiency*
Respirator Mask (PFE Standard)	USA: NIOSH (42 CFR 84)	N95 PFE ≥ 95% for 0.3-micron particles

* (PFE - Particle Filtration or BFE - Bacteria Filtration Efficiency)

Respiratory Mask of Class N95 is suitable for filtering of the most penetrating particles of 0.3-micron size that can be found in the air.


0.1 micron: Particle Filtration Efficiency standard (PFE).

0.3 micron: Used to represent the most-penetrating particle size (MPPS), which is the most difficult size particle to capture.

3.0 microns: Bacteria Filtration Efficiency standard (BFE) for masks per EN 14683.

		Boxes	Cases
ORDERING INFORMATION	Materials	350 gsm white cardboard	Brown 5-ply single wall
	Dimensions (in) (mm)	5.12*4.72*6.30 (in) 130*120*160 (mm)	24.41*11.02*13.38 (in) 620*280*340 (mm)
	8265	Gross Weight: (lb) 11.90lb/cs	(kg) 5.4kg/cs
Package 20 masks/polybag/box * 20 boxes = 400 masks / case			

ONE-CARE®

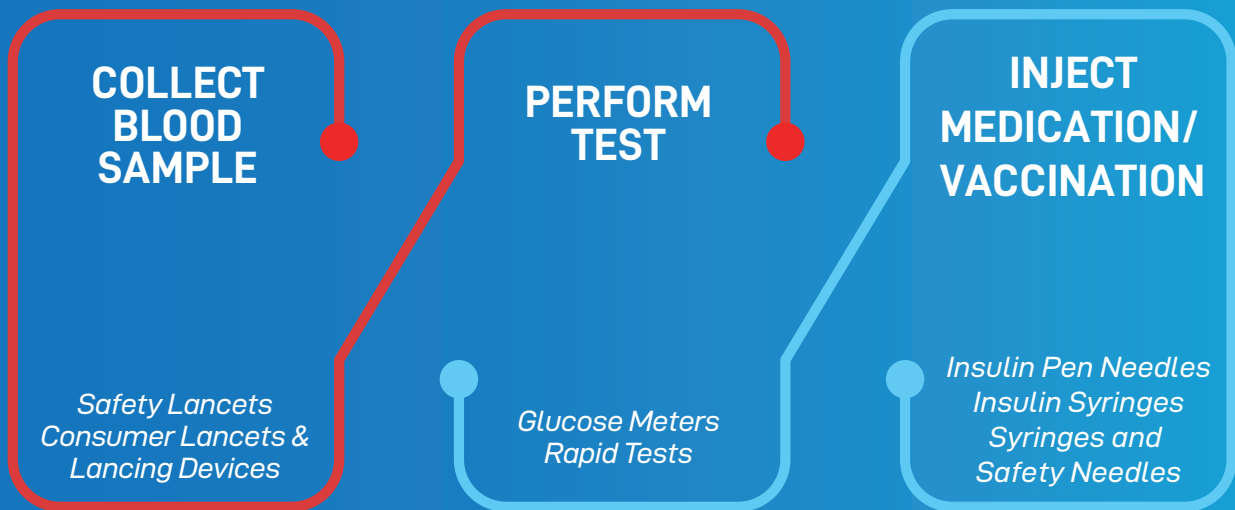

www.one-care.com

Americas

USA

📍 MediVena
1 Glenlake Parkway, Ste 700
Atlanta, GA 30328, USA
☎ +1-404-422-1513
✉ info@medivena.com
🌐 www.medivena.com

Europe

Poland

📍 MediVena Sp. z o.o.
Al. Jana Pawła II 27
00-867 Warsaw, Poland
☎ +48 661 920 142
✉ info@medivena.com
🌐 www.medivena.com