

ONE-CARE™

Lancing Device & Lancets



Blood micro-sampling has never been easier...

Easy
Easy Depth Setting
with clear marking

"Ready to Activate"
Visual Indicator

Compact Size
Compact Size for discrete
blood sampling

Convenient
Convenient Loading System



Universal
Universal Design for
fit with most common lancets



**Ultra-thin, electro-polished
lancets for comfortable sampling**

EASY

- Low fore to load the device with convenient end grip
- Adjustable depth setting easy to rotate and with clear visual marking

GENTLE

- Compact size with tight linear lancet guiding for gentle puncture
- Quite and soft activation for seamless blood sampling

SAFE

- Color-Coded loading indicator to distinguish when the device is ready to be used
- Automated ejector to safely remove lancet from the device

UNIVERSAL & ULTRA-FINO

- 23G to 33G ultra-thin lancets of universal fit
- Tri-Bevel, electro-polished lancets for enhanced user comfort




**Automated Ejector
for safe lancet removal**



INSTRUCTIONS FOR USE

- Intuitive and easy to follow steps for user comfort.

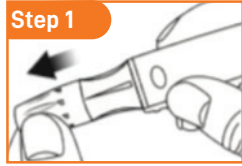


Remove the Device Cap




Insert a new Lancet Into the Device




Twist off and pull out safety cap of the Lancet

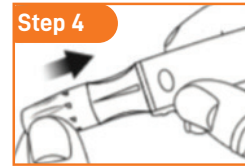


Put the Device Cap back on

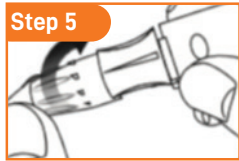

Adjust the lancing depth

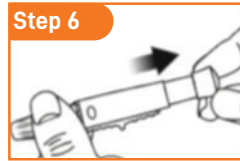

Load the Device (Click Indicator)

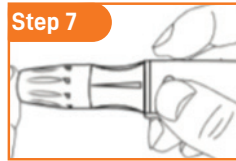

Activate the Device to obtain blood sample



Drive the Lancet into the safety cap



Eject the lancet safely from the device

ORDERING INFORMATION

CODE	DIAMETER (G) / DEPTH (mm/in)	DESCRIPTION	PACKAGING
1200	10 - Depth Setting	Universal Lancing Device	1 pc / box 100 pcs / case
1201	23G	23G Universal Blood Lancets	100 pcs / box 10.000 pcs / case
1202	26G	26G Universal Blood Lancets	100 pcs / box 10.000 pcs / case
1203	28G	28G Universal Blood Lancets	100 pcs / box 10.000 pcs / case
1204	30G	30G Universal Blood Lancets	100 pcs / box 10.000 pcs / case
1205	31G	23G Universal Blood Lancets	100 pcs / box 10.000 pcs / case
1206	32G	32G Universal Blood Lancets	100 pcs / box 10.000 pcs / case
1207	33G	33G Universal Blood Lancets	100 pcs / box 10.000 pcs / case

WWW.ONE-CARE.COM

