

ONE-CARE[®]

SAFETY PEN NEEDLES **OPTI**



Insulin delivery has never been easier...



Safety

Clear safety shield automatically locks after use to protect the needle

Simple

Color indicator becomes invisible after engagement of safety mechanism

Ergonomic Design

Ergonomic main housing with ribs for comfortable grip

Universal

Compatible fit with most pen injectors

SAFE

- Automated shielding mechanism ensuring that needle is covered before and after use
- Single-use mechanism preventing sharps injuries and reuse of the device
- All safety pen needles are ETO-sterilized and protected with cover and sterility peel-off tab

SIMPLE

- Safety mechanism engagement indicator via "Click" sound and disappearance of color coded inner shield after fully inserting the needle into the skin
- Optimally designed, clear shield for precise placement on injection site with improved visibility

GENTLE

- Design with concealed needle helps to relieve patient's anxiety
- Polished, silicone-coated, ultra-thin needles are ideally positioned during skin penetration thanks to precise tracks for the safety mechanism elements enhancing patients comfort
- Ribbed housing provides a comfortable grip for easy attachment and removal from the pen

UNIVERSAL

- Universal design for compatible use with most injection pens available on the market
- Tight fit prevents leakages of insulin
- Easy to recognize color-coded product versions address a variety of patient skin types and injection requirements

UNIVERSAL STANDARDS

- Rigorously tested per international standards for fit, functionality and reliable performance

Compatibility

FlexTouch® Novo Nordisk
 FlexPen® Novo Nordisk
 NovoPen® Echo Novo Nordisk
 NovoPen® 4 Novo Nordisk
 NovoPen® 5 Novo Nordisk
 Lantus SoloStar®
 OptiPen® Pro 1 Sanofi
 KwikPen® Humalog


Liprolog® Eli Lilly
 Huminsulin® Eli Lilly
 HumaPen® Luxura™ Eli Lilly
 HumaPen® Luxura™ HD Eli Lilly
 HumaPen® MEMOIR™ Eli Lilly
 Exenatide Pen Byetta
 Autopen® Classic Owen Mumford
 Autopen® 24 Owen Mumford

Omnican® Pen 31B. Braun
 Ypsopen® Ypsomed
 OptiClik® Ypsomed
 BerlinPen® Areo Berlin-Chemie
 BerlinPen® Areo 2 Berlin Chemie
 BerlinPen® 302 Berlin Chemie
 Tactipen® Haselmeier
 and many others

ONE-CARE™
 SAFETY PEN NEEDLES

INSTRUCTIONS FOR USE

- Intuitive handling with simple steps. Insulin delivery has never been easier...

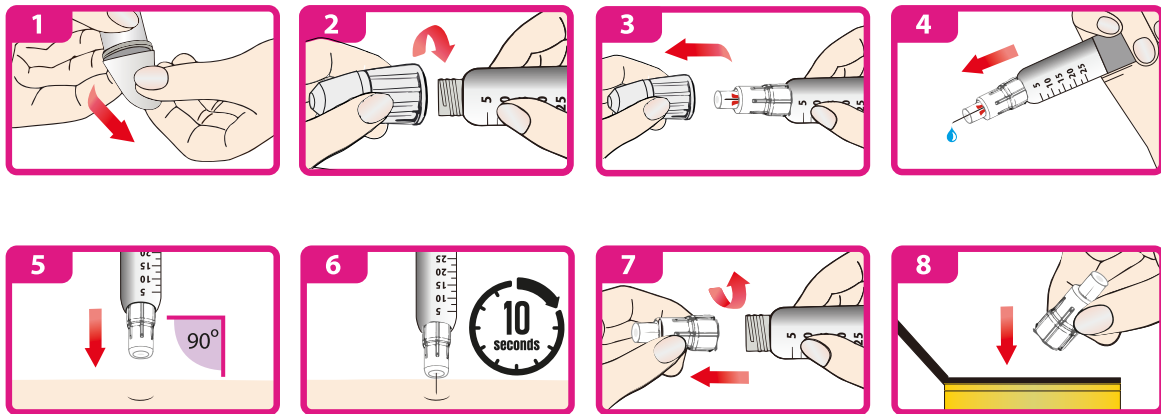

BEFORE
 USE



ONE-CARE® Pen Needle with
 sterility tab and outer cover

AFTER
 USE



ONE-CARE® Pen Needle after use
 (color-coded inner shield invisible)

PRODUCT SELECTION CHART

Code	Needle	Pcs per Box	Case Qty	Case Dimensions	Case Weight	Pallet Qty	Pallet Dimension	Pallet Gross Weight
3301	30G x 5mm	100	20 bx	13.6x12.6x9.1" 34.5x32x23 cm	16.31 lb 7.40 kg	72 cs 1,440 bx	47.3x39.4x78.8" 100x120x200 cm (2.4 CBM)	1,190.5 lb 540 kg
3302	30G x 8mm	100	20 bx	34.5x32x23 cm	7.40 kg	1,440 bx	47.3x39.4x78.8" 100x120x200 cm (2.4 CBM)	1,190.5 lb 540 kg

WWW.ONE-CARE.COM

MediVena
 inspired by innovation

Americas

USA
 MediVena
 1 Glenlake Parkway, Ste 700
 Atlanta, GA 30328, USA
 +1 404 422 1513
 info@medivena.com
 www.medivena.com

Europe

Poland
 MediVena Sp. z o.o.
 Al. Jana Pawla II 27
 00-867 Warsaw, Poland
 +48 22 299 44 19
 info@medivena.com
 www.medivena.com

