









ONE-CARE[™] PRO

GLUCOSE MONITOR & STRIPS
FOR PROFESSIONAL USE



Glucose monitoring has never been easier...

 FREE	 0.5µl	 AST	 14
Auto Calibration	0.5 µL Blood Sample	Alternative Site Testing	Instant 14 Days Average
 5 ^{sec}	 500		
5 Second Test	500 Memory Sets	PC Download	Strip Ejector



Bigger Buttons in a compact device

Easy to use



USB and Bluetooth (Optional)

Provide both wire and wireless connection for data transmission



Test Strip Ejector

Reduce the possibility of cross-contamination due to bloodborne pathogens



Test Strip

Longer shelf life 24 months



SPECIFICATIONS




Auto Coding

Simplified way to test, easy to use, no more misreading due to incompatible code card

GOD Enzyme


0.5 µl blood sample
Highly Accurate for Capillary Whole Blood samples

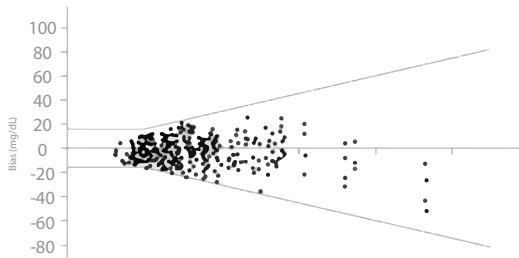

Hematocrit	30-55%
Test Sample	Capillary Whole Blood
Measuring Unit	mg/dL or mmol/L
Measuring Range	20 - 600 mg/dL (1,1 - 33,3 mmol/L)
Average	14/07/30 days
Reminder Alarm	4 sets
External Output	UART with mini USB socket (Optional) Bluetooth (Optional)
Power Supply	2 AAA batteries
Battery Life	Approximately 1,000 tests
Dimensions	100 (L) x 57 (L) x 24 (A) mm
Weight	58 g without batteries


INSTRUCTIONS OF USE

1 STEP

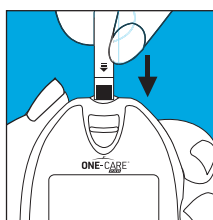


Insert a Test Strip

2 STEP

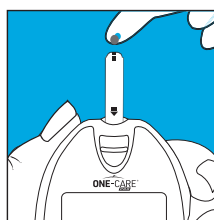


Apply Blood Sample

3 STEP

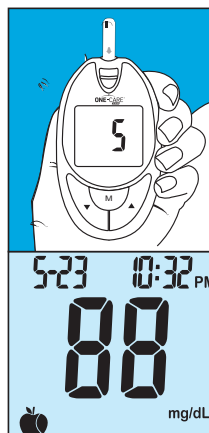

Wait 5 Seconds and Check Results

4 STEP

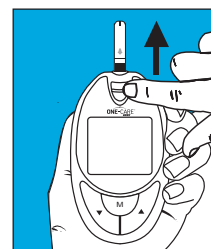


Discard the used Test Strip

ORDERING INFORMATION

CODE	PRODUCT	PCS per BOX	BOXES per CASE
4001	ONE-CARE™ PRO GLUCOSE METER FOR PROFESSIONAL USE	1	80
4002	ONE-CARE™ PRO GLUCOSE STRIPS FOR PROFESSIONAL USE	50	560
4003	ONE-CARE™ PRO CONTROL SOLUTION LEVEL 1	1	900
4004	ONE-CARE™ PRO CONTROL SOLUTION LEVEL 2	1	900
4005	ONE-CARE™ PRO CONTROL SOLUTION LEVEL 3	1	900
4010	ONE-CARE™ PRO CONTROL SOLUTION LEVEL 1 & 3 (Low & High)	2	900

WWW.ONE-CARE.COM

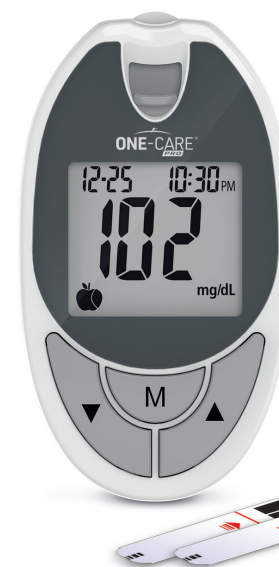

Americas

USA
MediVena
1 Glenlake Parkway, Ste 700
Atlanta, GA 30328, USA
☎ +1 404 422 1513
✉ info@medivena.com
🌐 www.medivena.com

Europe

Poland
MediVena Sp. z o.o.
Al. Jana Pawla II 27
00-867 Warsaw, Poland
☎ +48 22 299 44 19
✉ info@medivena.com
🌐 www.medivena.com

