

# ONE-CARE<sup>®</sup> PRO

GLUCOSE MONITOR & STRIPS  
FOR PROFESSIONAL USE



Glucose monitoring has never been easier...

<b>Auto Calibration</b>	<b>0.5 µL Blood Sample</b>	<b>Alternative Site Testing</b>	<b>Instant 14 Days Average</b>
<b>5 Second Test</b>	<b>500 Memory Sets</b>	<b>PC Download</b>	<b>Strip Ejector</b>



## Bigger Buttons in a compact device

Easy to use



## USB and Bluetooth (Optional)

Provide both wire and wireless connection for data transmission



## Test Strip Ejector

Reduce the possibility of cross-contamination due to bloodborne pathogens



## Test Strip

Longer shelf life 24 months



## SPECIFICATIONS




### Auto Coding

Simplified way to test, easy to use, no more misreading due to incompatible code card

### GOD Enzyme


0.5 µl blood sample  
Highly Accurate for Capillary Whole Blood samples

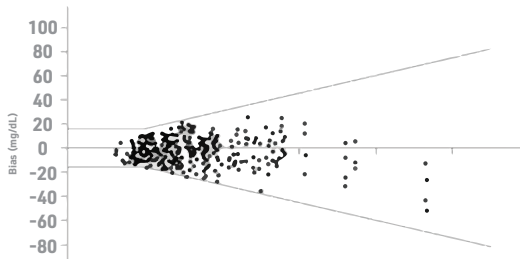

Hematocrit	30-55%
Test Sample	Capillary Whole Blood
Measuring Unit	mg/dL or mmol/L
Measuring Range	20 - 600 mg/dL (1,1 - 33,3 mmol/L)
Average	14/07/30 days
Reminder Alarm	4 sets
External Output	UART with mini USB socket (Optional)   Bluetooth (Optional)
Power Supply	2 AAA batteries
Battery Life	Approximately 1,000 tests
Dimensions	100 (L) x 57 (L) x 24 (A) mm
Weight	58 g without batteries


## INSTRUCTIONS OF USE

### 1 STEP

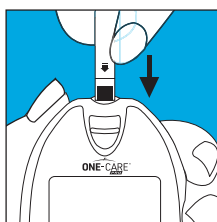


Insert a Test Strip

### 2 STEP

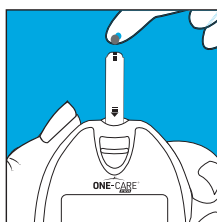


Apply Blood Sample

### 3 STEP

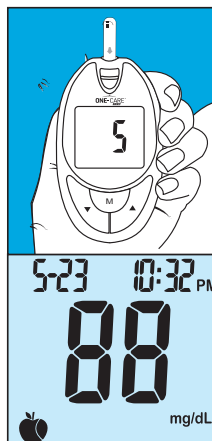

Wait 5 Seconds and Check Results

### 4 STEP

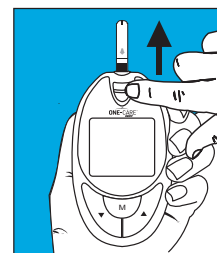

Discard the used Test Strip



## ORDERING INFORMATION

CODE	PRODUCT	PCS per BOX	BOXES per CASE
4001	ONE-CARE® PRO GLUCOSE METER FOR PROFESSIONAL USE	1	80
4002	ONE-CARE® PRO GLUCOSE STRIPS FOR PROFESSIONAL USE	50	560
4010	ONE-CARE® PRO CONTROL SOLUTION LEVEL 1 & 3 (Low & High)	2	900



[WWW.ONE-CARE.COM](http://WWW.ONE-CARE.COM)

#### Americas

**USA**  
MediVena  
1 Glenlake Parkway, Ste 700  
Atlanta, GA 30328, USA  
+1 404 422 1513  
info@medivena.com  
www.medivena.com

#### Europe

**Poland**  
MediVena Sp. z o.o.  
Al. Jana Pawla II 27  
00-867 Warsaw, Poland  
+48 22 299 44 19  
info@medivena.com  
www.medivena.com

